




CRANIAL PERFORATORS

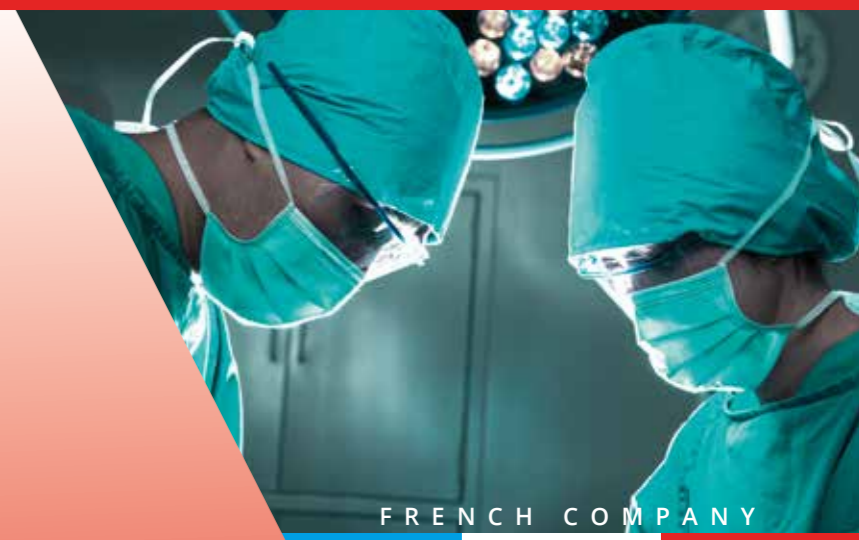

SINGLE-USE PRODUCT

FRENCH COMPANY

Vitalys Surgical



FRANCE & EXPORT HEAD OFFICE
193 Bd de Laval, 35500 VITRÉ - FRANCE
Phone: +33 (0)2 90 64 00 40
Fax: +33 (0)2 90 64 00 42

EXPORT SALES DEPARTMENT
Phone: +33 (0)2 90 64 00 41

FRANCE SALES DEPARTMENT
Phone: +33 (0)2 90 64 00 40

Ref. 0106-00-00 - 01/01/2020
© grenandwhite.fr - Printed in France.



ADULT



ADOLESCENT



PAEDIATRIC

CE 0459  ASKORN MEDICAL - Phone: +33 (0)2 23 35 53 35 - 44c rue de Bray, 35510 CESSON-SÉVIGNÉ - FRANCE

DISTRIBUTED BY

www.vitalys-surgical.com

www.vitalys-surgical.com

Vitalys Surgical is a French company based in Brittany specialising in neurosurgery.

All the products developed by the company have been **researched and tested by engineers** from a surgical-sector ISO 13485-approved specialist design consultancy.



Vitalys Surgical proposes a range of **100% MADE IN FRANCE** cranial perforators, available in **3 sizes to fit each morphology:**

ADULT

ADOLESCENT

PAEDIATRIC

> Serving neurosurgeons

Vitalys Surgical's founder, Pierre-Yves Collerai, has built up extensive experience working alongside neurosurgeons in the operating theatre and in the implantable medical devices sector. He has developed the company's products based on his observation of surgeon's practices. Its range of innovative cranial perforators is designed to **optimise the practitioner's procedures, comfort and safety.**

> Cutting-edge R&D

The medical devices proposed have been specially designed by R&D to provide neurosurgeons with a **completely safe, convenient-to-use range of products suitable for all cranial anatomies.** A **patent** has been filed to protect the **engage/disengage** concept of the VitalDrill cranial perforator.

> A safe secure process

The perforator's safety is primarily based on:

- The **anti-skid centring tip** at the start of perforation;
- The **optimised milling** of the perforator's cutting teeth;
- **Immediate release** at the right moment at the end of each perforation;
- The **formation of a bony capsule** at the end of perforation, creating a protective shield between the perforator drill and the dura mater. This safety guarantee is practical and reassuring for practitioners, who thus retain total and systematic control of their surgical procedure.

Each perforator is **controlled 5 times**, with activation of its mechanism, before being packaged in a blister.

> France & International Department

Vitalys Surgical and its sales team are at your disposal to handle all your requests, **in France or internationally.**

All our cranial perforators have an **auto-release mechanism** and are **single-use.**





SINGLE-USE CRANIAL PERFORATORS

ADULT SIZE Diameter 14/11

Reference 0106-03-00




VitalDrill ADULT

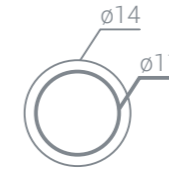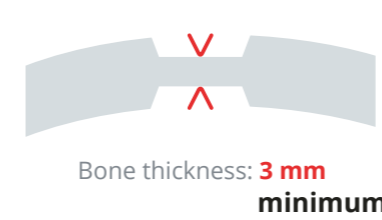

ADOLESCENT SIZE Diameter 11/7

Reference 0106-02-00

VitalDrill ADO

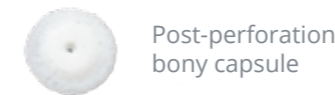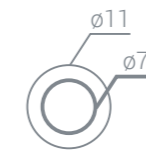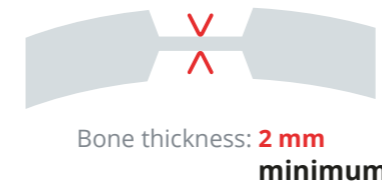

PAEDIATRIC SIZE Diameter 9/6

Reference 0106-01-00

VitalDrill CHILD

